

PLUNGER TYPES




Fitted Plunger

Ideal syringe used for dispensing small volume liquids. Plunger is individually fitted to glass barrel and cannot be replaced.




Gastight (GT) - PTFE

Plunger has a pure PTFE tip creating a strong seal inside the syringe barrel for dispensing of liquids and gases. PTFE tip has good solvent compatibility. Plunger is replaceable.

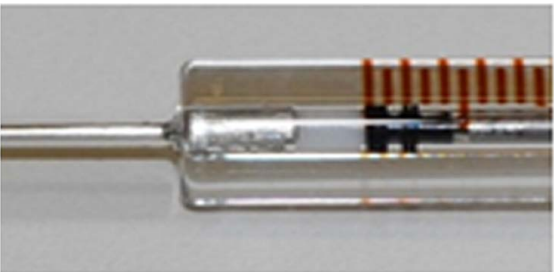

Gastight (GT) - PTFE - Black Seal

Plunger has a Black hybrid PTFE tip creating a strong seal inside the syringe barrel for dispensing of liquids and gases. The hybrid PTFE seal has broader solvent resistance. Plunger is replaceable.

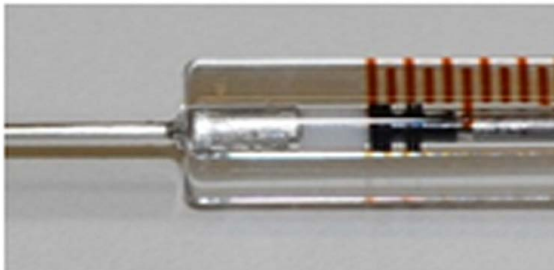

Gastight (GT) - PE Seal

Plunger has a Polyethylene (PE) tip creating a strong seal inside the syringe barrel for dispensing liquids and gases. PE is good for buffers and aqueous solutions. Plunger is replaceable.